


## Pillendose Diagramm



## Vor Gebrauch!

- We recommend washing each pill box before using your Medinizer for the first time. You can wash them by hand with a mild detergent or in a dishwasher at 50 °C.
- We also recommend cleaning before refilling your pill boxes, especially if your medication changes
- Remember to keep your pill boxes out of reach of from children.
- Pay close attention to the storage requirements of your medication; some medication cannot be removed from its original packaging before use. There may be temperature or photo-sensitivity or other requirements. If in doubt please ask your doctor or pharmacist

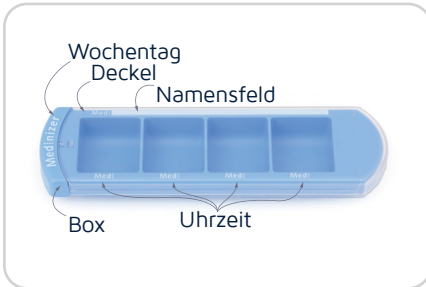

- Remove the pill boxes from the cover
- Slide the lid/lids back and remove completely



- Place your medication in the appropriate compartments/boxes



- Replace the lid. The lid is secure when you feel/hear a click



- Schieben Sie den Deckel zurück, sodass er mit dem entsprechenden Teiler ausgerichtet ist
- Geben Sie das Medikament einfach in die Hand



- Schieben Sie den Deckel in die geschlossene Position
- Trinkgeld!** Auf der mitgelieferten Medikamentenkarte können Sie Informationen zu Ihren Medikamenten vermerken



- Trinkgeld!** Unter dem Boden jedes Fachs befinden sich 1-4 hervorstehende Punkte für die verbesserte Nutzbarkeit von sehbehinderten Menschen